

# Hoist Accessories



## Jordan Frame (& Strap Kit)

- Supine lifting device including 2 piece frame, gliders & strap kit for hoist
- 2 sizes – wide (SWL 200kg) and extra wide (SWL 230kg)
- Suitable for yoke hook attachment (preferably 4 Point Yoke)

**Available on:** Alto, Kratos and other large non-Allegro hoists

**Indications:** Lifting patients requiring spinal or orthopaedic immobility

**Useful for:** Floor lifts where spinal or orthopaedic injury is suspected



Jordan Frame - LF01

## Weigh Scale

- Universal weighing module with tare function and kg/lb conversion
- SWL up to 300kg
- Suitable for yoke or pivot attachments

**Available on:** All Allegro & Allegro Lite hoists plus some non-Allegro hoists

**Indications:** Where weighing of patients is difficult due to the level of support required

**Useful for:** Patients requiring constant monitoring of weight for general wellbeing



Weigh Scale - WD03

**Allegro**  
CONCEPTS PTY LTD

Innovation in Patient Handling

## Hoist Accessories

Allegro Concepts offer a broad range of accessories that interface with the Allegro steel frame lifting hoists, the Allegro Lite aluminium frame lifting hoists plus many non-Allegro lifting hoists.

Some accessories (such as the Jordan Frame) are designed to increase the functionality of existing patient lifting equipment. In this case a regular patient lifting hoist with yoke hook attachment can be used to lift patients requiring total spinal immobility e.g. unstable quadriplegia immediately post trauma.

### Sling Bag

- Heavy duty vinyl bag for storing slings and other patient handling accessories (e.g. Swift Slide Sheets)
- Attaches to the push handles of the hoist with 2 Velcro tabs
- Custom sizes and monogramming available on request

**Available on:** All Allegro and non-Allegro patient lifting hoists

**Indications:** Where storage of and immediate access to the sling is required

**Useful for:** Storing slings, Swift Slide Sheets and Patient Transfer Belts on the end of patients beds



Sling Bag - SB01

### Tempo Handle Kit

- Pair of bolt on handles with non slip rubber grip
- Designed to provide optimum grip position on Salsa standing transfer hoist
- Multi adjustable in angle, height and distance from patient

**Available on:** Allegro Salsa standing transfer hoist

**Indications:** Where the regular grip position is not comfortable or ergonomically placed for the patient

**Useful for:** Hemiplegic patients who require optimum gripping position with the one functional hand



Tempo Handle Kit - THK01

### Leg Strap Knob Kit (Salsa 200)

- Hardware for the attachment of leg straps when custom fitted to an Allegro standing sling
- Fitted to upper clevis of lifting actuator (can be retrofitted)
- Allows the use of the Allegro Salsa for non-standing transfers

**Available on:** Allegro Salsa standing transfer hoist

**Indications:** Patients with sitting balance but limited (or fluctuating) weight-bearing capacity

**Useful for:** Conditions such as Multiple Sclerosis where transfer ability can fluctuate within a short space of time



Knob Kit - KK01


### Leg Protector

- Clear Polycarbonate plastic covers for protecting paint finish
- Fitted and removed without tools
- Easily cleaned for infection control

**Available on:** Allegro steel frame patient lifting hoists

**Indications:** Facilities which rely heavily on patient lifting hoists and are concerned with general upkeep of equipment

**Useful for:** Preventing rusting of patient lifting hoists used in wet environments

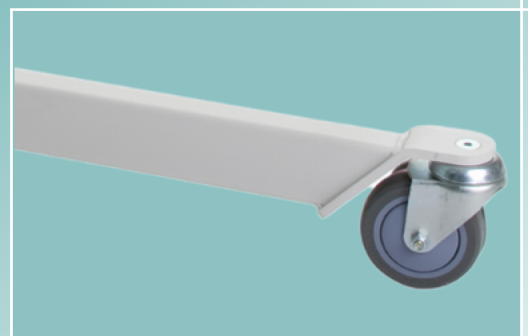

Leg Protector - LPA01

## Allegro Concepts Pty Ltd

Unit 9, 110-120 Silverwater Road | Silverwater | NSW | 2128

☎ 02 6895 2888 | ✉ sales@allegroconcepts.com.au

🌐 www.allegroconcepts.com.au