



www.jbmedical.com.au

JANUARY 2023

BROCHURE


WIRELESS FLOOR SENSOR MAT

FLOOR MAT & TRANSMITTER: MED-FM-03X-WT

FULL KIT: RECEIVER BOX, FLOOR MAT & TRANSMITTER: MED-FM-03K-WT(X)

Indigo Care provides a smart Australian designed 2-Fold Premium WIRELESS Floor Mat. The cordless floor mat is purposely designed for easier storage and trip-hazard free. The transmitter is an exclusive ID encoded RF Transmitter which communicates with the Wireless Receiver Box (RRU) or the Universal Call Point (UCPV5.1). It allows the Floor mat to wirelessly transmit an alarm to the call point. Simply connect the Transmitter to the floor sensor mat. Program the Transmitter to RRU or Dual Function Call Point by following the instructions printed. The sensor mat will now trigger the call point when activated within the transmitting area. (Note: Transmitter and Floor Mat can be ordered separately).

12 MONTHS WARRANTY

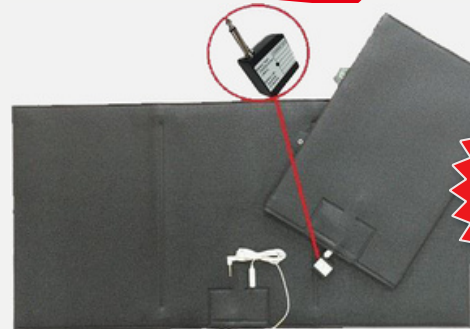

\$ XXXX
ex gst

KEY SPECIFICATIONS:

Floor mat

- Dimensions: 610mm X 1220mm X 6mm (2-Fold)
- Electrical Resistance < 100 Ohm
- Black PVC Yoga Pattern
- Non-slip bottom surface
- Operation Weight: 30-200kg

Transmitter

- Connection: 6.3mm Jack socket
- No External Antenna
- LED low battery indicator
- Transmits up to 30m (line of sight)
- Frequency: 433.92MHz
- Uses a 3V battery-CR2032 which will last up to 12 months
- Warranty: 12 months.

KEY FEATURES

- ✓ 2-FOLD MAT DESIGN FOR SAFETY, EASY STORAGE AND MULTIFUNCTION
- ✓ REDUCE THE RISK OF FALLS & INJURIES
- ✓ NO WIRES- NO TRIP HAZARDS
- ✓ LOW PROFILE, NON-SLIP MATERIAL
- ✓ WATERPROOF & WASHABLE
- ✓ PORTABLE
- ✓ EASY TO USE
- ✓ EASY TO MOVE FROM ROOM TO ROOM
- ✓ INTEGRATES TO MOST NURSE CALL SYSTEMS
- ✓ PLUG & PLAY OPERATION
- ✓ DURABLE & RELIABLE
- ✓ COST EFFECTIVE



ON SALE
ORDER NOW

WIRELESS CUSHION FLOOR SENSOR MAT

CUSHION FLOOR MAT & TRANSMITTER: MED-FM-04X-WT

FULL KIT: RECEIVER BOX, FLOOR MAT & TRANSMITTER: MED-FM-04K-WT

Indigo Care provides a smart Australian designed Foldable Premium WIRELESS Cushion Floor Mat for residents who are at risk of falling out of their bed. The cordless floor mat purposely designed for easier storage and trip-hazard free. The transmitter is an exclusive ID encoded RF Transmitter which communicates with the Wireless Receiver Box (RRU) or the Universal Call Point (UCPV5.1). It allows the Floor mat to wirelessly transmit an alarm to the call point. Simply connect the Transmitter to the floor sensor mat. Then teach in the Transmitter to RRU or Dual Function Call Point by following the instructions printed. The sensor mat will now trigger the call point when activated within the transmitting area. (Note: Transmitter and Floor Mat can be ordered separately)



\$ XXXX
ex gst

KEY FEATURES

- ✓ PREMIUM OUT OF BED CUSHION MAT
- ✓ THICK FIRM FOAM TO REDUCE, IF NOT ELIMINATE INJURIES
- ✓ INCLUDE A 2-PIECE SENSOR MAT AND A WIRELESS TRANSMITTER
- ✓ FOLDABLE FOR EASY STORAGE
- ✓ EASY TO MOVE FROM ROOM TO ROOM
- ✓ SPRAY AND WIPE TO CLEAN WITH ANTIBACTERIA SOLVENT
- ✓ INTEGRATES TO MOST NURSE CALL SYSTEMS
- ✓ PLUG & PLAY OPERATION
- ✓ COST EFFECTIVE
- ✓ LARGE PRESSURE MAT (1800MM X 750MM X 40MM)

KEY SPECIFICATIONS

- Dimensions: 1800mm X 750mm X 40mm (Foldable)
- Electrical Resistance < 150 Ohm
- Blue PVC Soft Foam Design
- 40 mm of Foam Padding
- Connection: 6.3mm Jack socket
- No External Antenna
- LED low battery indicator
- Transmits up to 30m (line of sight) Frequency: 433.92MHz
- Uses a 3V battery-CR2032
- Warranty: 12 months

PRODUCT ACCESSORIES

- Wireless Receiver (RRU) Wireless Bed Mat
- Wireless Under Mattress Mat Wireless Chair Mat
- Wireless Waterproof Neck or Wrist Pendant
- Wireless PIR Motion Sensor



NDIS REGISTERED SERVICE PROVIDER



1300 788 182